


Tecan SPECTRAFluor

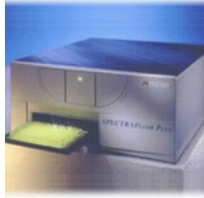

Manufacturer: Tecan
Model Number: SPECTRAFluor
Web Address: www.tecan.com

SPECTRAFluor measures any kind of plate formats (from Petri dish to 1536-well plates). It contains three measurement modes: fluorescence, absorbance and luminescence.

Plate Reader: Commands

- ▶ **GetServerVersion()** - Retrieve the version and licence information for the Magellan server.
- ▶ **PlateIn()** - Move the plate carrier into the instrument.
- ▶ **PlateOut()** - Move the plate carrier out of the instrument.
- ▶ **GetPlatePosition()** - Retrieve the current plate position.
- ▶ **SetFileName(fileName)** - Set the name of the output file for the method.

fileName	String	Name of the output file for the method.
----------	--------	---

- ▶ **SetEmpty(state)** - Reset the state of the workspace.

state	Boolean	True resets the workcell to the empty state.
-------	---------	--

- ▶ **InsertMethod(method)** - Load the specified method into the current workspace.

method	String	Name of the method to be loaded.
--------	--------	----------------------------------

- ▶ **Run()** - Run the previously loaded method.

Plate Reader: Errors

- 🔥 **SPECTRAFluorError(description)** - Error occurred during command execution.

description	String	Error description
-------------	--------	-------------------

ReTiSoft Inc.
 366 Revus Avenue, Unit 21
 Mississauga, Ontario, Canada, L5G-4S5
 Main: 647-724-2398 Europe: 33-9-7518-0225

Web: www.retisoft.ca Email: prodziew@retisoft.ca